# **Data & Instructions**


## Hydrazide-activated SepFast MAG

Hydrazide activated SepFast MAG is a magnetic agarose resin that is particularly useful for making glycoprotein-attached magnetic adsorbents.

Hydrazide-activated agarose permits the coupling of aldehyde- or ketone-containing ligands through the formation of stable hydrazone linkages. It is a particularly powerful way to immobilize glycoproteins via their carbohydrate chains while leaving critical active sites intact. It is a method called "site-directed immobilization". For example, antibodies can be immobilized to an agarose solid support through a linkage to their glycosylated side chain that leaves their Fab region free for antigen binding.

Glycoproteins can be readily oxidized with sodium periodate to generate formyl groups on their carbohydrate chains before the coupling reaction with Hydrazide-activated SepFast MAG. It is a successful and well-documented technique. The coupling reaction is mild and easy to carry out. No toxic chemical or special equipment is required.

Hydrazide-activated SepFast MAG resins can be readily employed to make custom affinity chromatography media for both small scale and large scale purification applications.

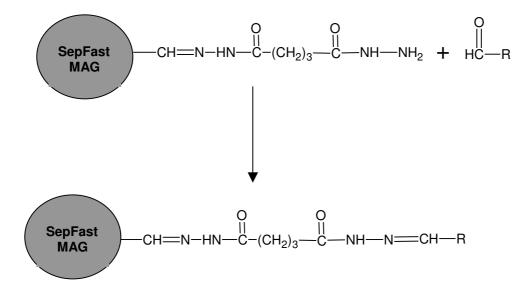

## 1. Properties

Hydrazide activated SepFast MAG 4HF is made of highly cross-linked 4% agarose containing magnetic material. It shows high mechanical rigidity and high chemical stability for harsher operational conditions.

Agarose has long been used for chromatographic separations due to its excellent hydrophilic and low non-specific-binding nature. The particles have an open pore structure with excellent mass transfer properties to large protein molecules.

The base matrix is activated with hydrazide through a long hydrophilic spacer arm. Hydrazideactivated SepFast MAG media is supplied as an aqueous suspension. The main characteristics are summarized in Table 1.

Group to be coupled	-HC=O (aldehyde), -C=O (ketone)	
Matrix	SepFast MAG 4HF: highly cross-linked 4% magnetic agarose	
Particle size*	50 – 150 μm	
Coupling level	>2 mg hlgG / ml medium	
pH stability	4 -11 (ligand dependent)	
Chemical stability	Compatible with all commonly used aqueous chemicals, provided the ligand to be coupled can withstand	
Storage	+4 <sup>o</sup> C - +8 <sup>o</sup> C	

## Table 1: Characteristics of Hydrazide-activated SepFast MAG 4HF:

\*Other particle sizes available on request.

#### 2. Ligand immobilization

#### 2.1 Oxidation of glycoproteins

Glycoproteins need to be oxidized with sodium periodate to generate formyl groups on their carbohydrate chains first.

2.1.1 Dissolve or dilute glycoprotein to a concentration of 1 - 10 mg / ml in coupling buffer, 0.1 M sodium phosphate, pH 7.0, in an amber container.

2.1.2 Add sodium meta-periodate (5 mg / ml coupling buffer). Gently mix to dissolve the oxidizing agent.

2.1.3 Incubate the sample in the dark at room temperature for 30 minutes. The reaction is light sensitive.

2.1.4 Stop the reaction by desalting and buffer exchange in a suitable gel filtration column.

#### 2.2 Coupling reaction

The following is a general ligand coupling procedure.

2.2.1 Wash Hydrazide-activated SepFast MAG medium with at least 5 gel volumes of coupling buffer in a filtration device.

2.2.2 Transfer the washed and suction dried gel to the oxidized glycoprotein solution. The volume ratio between protein solution and gel should be 1 : 1.

2.2.3 Mix the slurry at room temperature overnight.

2.2.4 Wash the gel with 5 volumes of 0.1 M Tris/HCl + 0.5 M NaCl, pH 8.0, followed with 5 volumes of 0.1 M acetate buffer + 0.5 M NaCl, pH 4.0. Repeat this washing cycle two more times.

2.2.5 Wash the gel with working / equilibration buffer before use.

## 3. General considerations over the immobilization efficiency

## 3.1 pH

A buffer with neutral pH is most frequently used for protein immobilization. However, the coupling pH may be optimized between 6 to 10 to get the best result (e.g. high coupling yield with high ligand activity).

## 3.2 Coupling solution

A solution containing ketone groups should be avoided.

Certain organic solvents in diluted format may be introduced to improve the solubility of the ligand. The suitability of such solvents should be tested in advance.

#### 3.3 Further reductive amination

The stability of the hydrazone bond formed between hydrazide and formyl groups can be further improved by mild reductive amination. Sodium cyanoborohydride can be added to convert the hydrazone bond to a secondary amino bond.

## 3.4 Washing of the final medium

The non-attached or weakly attached ligand needs to be fully washed away after the coupling reactions. A washing method employing alternating high pH and low pH can ensure an efficient removal of the unwanted species.

#### 4. Use of the immobilized affinity medium

The ligand coupled medium can be used for purifications using batch stirred tank mode or packed column mode. Handling of this material follows the same principles as handling of other agarose-based media.

## 5. Storage

Hydrazide-activated SepFast MAG media should be stored under 8°C. The coupled wet medium should be stored in the presence of a bacteria-proof agent (e.g. 20% ethanol) at 4-8°C. Never freeze the coupled medium.

#### 6. Further information

Visit www.biotoolomics.com for further information or contact the technical team or sales representatives.

#### 7. Ordering information

Product	Quantity	Code no.
Hydrazide-activated SepFast MAG 4HF	5 ml	450204
	50 ml	450205
	1 litre	450206



BioToolomics Ltd Unit 30A, Number 1 Industrial Estate Consett County Durham, DH8 6TJ United Kingdom

#### www.biotoolomics.com

Registered or registration-pending trademark of BioToolomics Litd: BioToolomics, SepFast. SuperSpin

All goods and services are sold subject to the terms and conditions of sale of BioToolomics Ltd. The company reserves the rights, subject to regulatory or contractual approval, if required, to make changes in the specifications and features shown herein, or discontinue the products described at any time without notice or obligation. Contact BioToolomics Ltd for the most current information.

@ 2010-2015 BioToolomics Ltd – All rights reserved.